


PIG® Lab & Chemical Cleanup Kit

KIT5007 Absorbs Acids; Bases; Unknowns, 1 each

Be prepared for small spills in the laboratory with our grab-and-go clean up kit for acids, bases and chemicals.

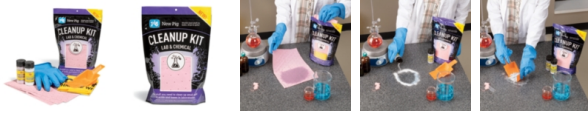

- Compact design lets you store a spill kit virtually anywhere you need one for fast response
- Easy-open, resealable bag protects kit contents from dirt and moisture
- Single-use kit is faster and more efficient than dragging a large spill kit to a small cleanup; safer and more effective than using rags or paper towels
- Includes hazmat absorbent pads, acid neutralizing loose absorbent and base neutralizing loose absorbent to clean up most acids, bases and chemicals (absorbs up to 25 ounces)
- Disposable gloves, mini dustpan and disposal bag also included to aid cleanup
- Only PIG Spill Kits feature PIG Absorbents - the world's #1 selling brand
- Spill Kit dimensions are approximate



Specifications

Container Type	Single Use Bag
Fluid Absorbed	Acids; Bases; Unknowns
Color	Purple
Dimensions	4" W x 9.25" L x 13" H
Absorbency	Up to 25.2 oz.
Brand	PIG
Sold as	1 each
Weight	1 lb.
# per Pallet	40
Composition	Mats - Polypropylene, Acid Neutralizer - Polymer, Sodium Carbonate, Alizarin, Base Neutralizer - Polymer, Citric Acid, Thymol Blue
Includes	3 - 10" W x 13" L PIG® HazMat Chemical Absorbent Mat Pad 1 - 3.5 oz. Acid Encapsulating/Neutralizing Absorbent 1 - 3.5 oz. Base Encapsulating/Neutralizing Absorbent 1 - Disposable Scoop and Scraper

1 - 9.25" X 13" x 4" Bag Single Use - Lab & Chemical Cleanup Kit
1 - 12" L SHOWA NITRI-DEX® Nitrile Gloves (GLV970-XL)
1 - 18" W x 30" H Polyethylene Disposal Bags (BAG201-S)

UNSPSC 47131905

UV Resistant No

Pigalog® Page Number Page 120

Metric Equivalent

Absorbency Up to 745.2 mL

Dimensions 10.2cm W x 23.5cm L x 33cm H

Weight .5 kg

Technical Information

CountryNoShip

This product is subject to export restrictions and cannot be shipped to: Canada, European Union, Great Britain.

Technical Documents

[PIG Hazardous Chemical Absorbents](#)

[Acid Neutralizer & Absorber](#)

[Base Neutralizer & Absorber](#)

[PIG® HAZ-MAT Absorbents](#)

Disclaimers

Acid Absorbents Notice

Not for use with hydrofluoric acid as it may create a hazardous gas.

Base Absorbents Notice

Avoid contact with metallic nitrates, cyanides, sulfides and strong oxidizers. Contact with sodium or calcium hypochlorite creates chlorine gas.

HAZ-MAT Notice

To ensure effectiveness and your safety, we recommend that you conduct compatibility and absorption testing of your chemicals with PIG® HAZ-MAT products prior to purchase. If you have any questions or need samples to test, please call Technical Services.



One Pork Avenue • Tipton, PA 16684-0304

1-855-493-4647 • Fax: 1-800-621-7447 • newpig.com • hothogs@newpig.com

